



**EZ-AF9000 - shown with E-Z Heat Controller and Heated Surgery Beds**

### **EZ-AF9000 Auto Flow System**

Precise, consistent, preset flow rate from each outlet, no variables...

The E-Z Anesthesia Auto Flow System offers highly efficient, exact gas flow that is easy to use.

- Ideal for work with multiple animals and unique configurations
- Precise, low flow gas saves money and protects personnel
- Flow rates are set specifically for mice and rats
- Five outputs supply anesthesia to the induction chamber and up to four breathers
- Preset flow rates eliminate the need for flow meters and manual adjustments
- Protocols are easily followed without calculations
- Flow rate is not impacted by the type or the number of accessories connected
- Minimal opportunity for operator error
- Each output is individually controlled with a simple on/off switch
- Pushbutton oxygen flush purges induction chamber before removing animals
- Carry handle facilitates safe, easy transport of the unit

### **EZ-AF9000-A Customized Auto Flow System**

- Outputs can be individually customized to specific flow rates
- Can be customized to accommodate guinea pigs and rabbits

#### **Specifications:**

**Anesthetic: Isoflurane O - 5%**

**Carrier Gas: Oxygen or Air**

**Flow Rate: Chamber - 1 LPM each**

**Breathers - 0.5 LPM each**

**Number of Outputs: 5**

Your Distributor is:  
Braintree Scientific, Inc.  
PO Box 850498, Braintree, MA 02185  
781-917-9526  
Email: [Info@braintreesci.com](mailto:Info@braintreesci.com)  
Web: [www.braintreesci.com](http://www.braintreesci.com)

## EZ-AF9000 Base System Includes:

EZ-190F Auto-Flow Anesthesia Machine

(1) Oxygen or Air Regulator

(2) EZ-103A Microflex Breathers

(1) EZ-178 Sure-Seal Mouse/Rat Induction Chamber


(1) Case of ReFresh Charcoal Filters (8)

Also includes mouse/rat nosecones, nosecone plugs, plastic surgery bed, and all connecting tubes

## EZ-AF9000 Base System



### EZ-AF9000 Auto Flow System

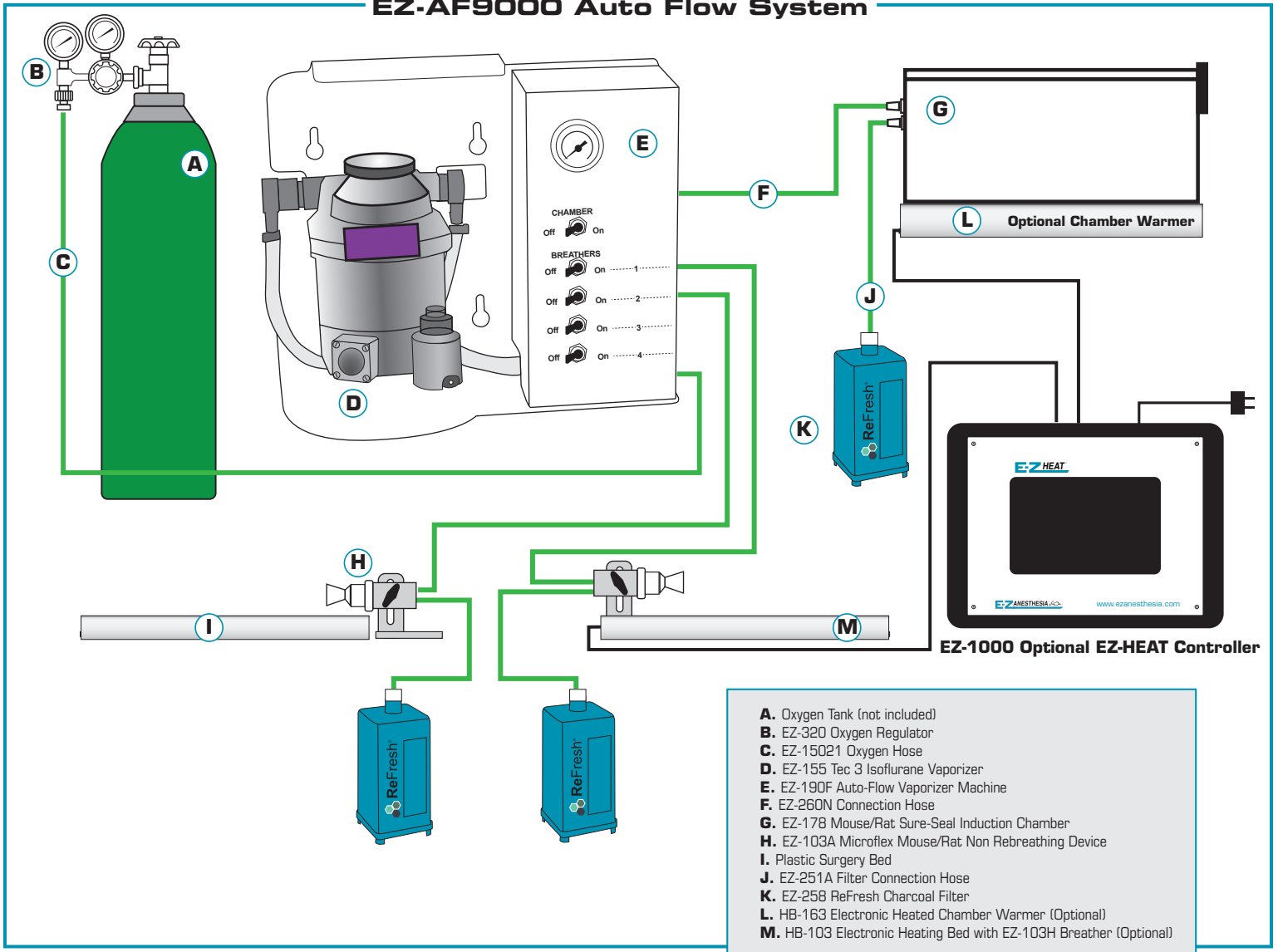

- A. Oxygen Tank (not included)
- B. EZ-320 Oxygen Regulator
- C. EZ-15021 Oxygen Hose
- D. EZ-155 Tec 3 Isoflurane Vaporizer
- E. EZ-190F Auto-Flow Vaporizer Machine
- F. EZ-260N Connection Hose
- G. EZ-178 Mouse/Rat Sure-Seal Induction Chamber
- H. EZ-103A Microflex Mouse/Rat Non Rebreathing Device
- I. Plastic Surgery Bed
- J. EZ-251A Filter Connection Hose
- K. EZ-258 ReFresh Charcoal Filter
- L. HB-163 Electronic Heated Chamber Warmer (Optional)
- M. HB-103 Electronic Heating Bed with EZ-103H Breather (Optional)

**BRAINTREE**  
**SCIENTIFIC, INC.**  
 LAB RESEARCH PRODUCTS

Your Distributor is:  
 Braintree Scientific, Inc.  
 PO Box 850498, Braintree, MA 02185  
 781-917-9526  
 Email: [Info@braintreesci.com](mailto:Info@braintreesci.com)  
 Web: [www.braintreesci.com](http://www.braintreesci.com)